

P/O NO.:		Inspector & Date	Shi Miao	15 th July. 2020	
Lot No.:	GL200610				
Supplier branch					
Supplier:	Green	Inspection Type:	Pre-shipment		
Inspection method:	Level S-4 for Workmanship appearance, S-1 for specific inspection, AQL 4.0				
Comment:	1. Based on the sampling, over quality is acceptable				

Item No.	Release Quantity	Carton Quantity	Finished Quantity	Unit in Pack	Unit in Carton	GW	Size of carton
							cm
Vinyl Glove-M	1000,000 pcs	1000 ctns	70 %	100 pcs	10 boxes	4.8kg	31.5*23*23
Vinyl Glove-L	1200,000 pcs	1200 ctns	80 %	100 pcs	10 boxes	5.3kg	31.5*23*23
Vinyl Glove-XL	1200,000 pcs	1200 ctns	60 %	100 pcs	10 boxes	5.8kg	31.5*23*23

Item no.	Sampling Size	Sampling Carton Qty.	Allowed Defects			Founding Defects			Specific Test
			Cr.	Ma.	Mi.	Cr.	Ma.	Mi.	
Vinyl Glove	125 pcs	12 ctns	0	7	10	0	0	0	5 pcs/size

Checklist	Checked	N.C.	N/A.	Remark/Comment/Result
Condition/Shipping mark/code	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Barcode/Sticker label/Packing	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Product style/Design/Artwork	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Workmanship/Appearance	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Dimension/Specification	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Country of origin	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Warning Instruction	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Safety requirement	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Physical test	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	


Workmanship & Appearance Defective Finding				
Defect description	Critical	Major	Minor	Defect photo no.
Total defectives				/
Maximum allowed	0	7	10	/




Dimension & Functional test							
Dimension-M	Spec	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Result
Total length (mm)	230+-5	235	235	235	235	235	Pass
Palm width (mm)	95+-5	95	95	95	95	95	Pass
Weight (g/pc)	4.0+-0.2	4.08 g(average of 10pcs)					Pass
Dimension-L	Spec	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Result
Total length (mm)	230+-5	235	235	235	235	235	Pass
Palm width (mm)	105+-5	105	105	105	105	105	Pass
Weight (g/pc)	4.5+-0.2	4.48 g(average of 10pcs)					Pass
Dimension-XL	Spec	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Result

Total length (mm)	230+5	235	235	235	235	235	Pass
Palm width (mm)	115+5	115	115	115	115	115	Pass
Weight (g/pc)	5.0+0.2	4.92 g(average of 10pcs)					Pass
Functional Test	Spec	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5	Result
Wearing comfort	Stimulation	ok	ok	ok	ok	ok	Pass
Open smoothly	In 1 second	ok	ok	ok	ok	ok	Pass
Elongation (%)	320	352	362	373	370	362	Pass
Tensile Strength (MPa)	9	11.21	12.35	10.25	10.67	14.22	Pass
Force at break (N)	/	3.8	4.6	4.2	4.6	4.5	Pass
Water leakage	No leaks	No leaks	No leaks	No leaks	No leaks	No leaks	Pass

Product Photos



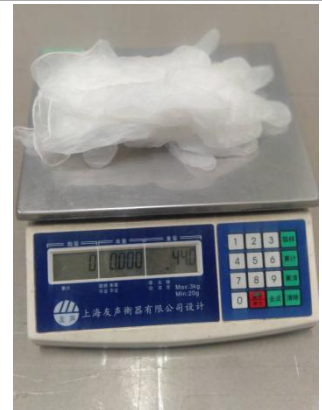

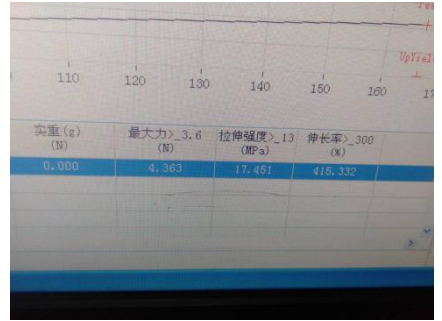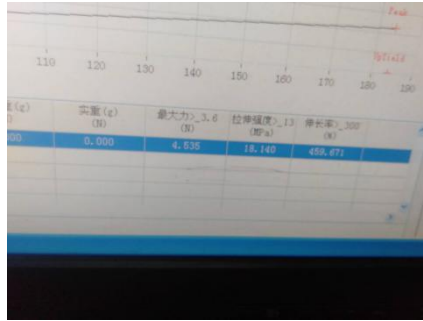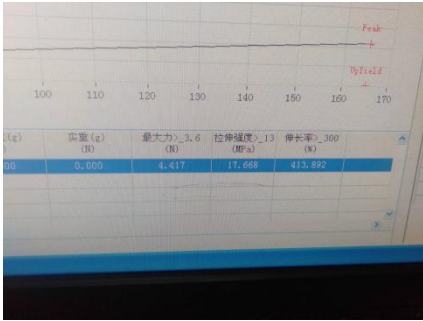

Defect Photos

--	--	--