
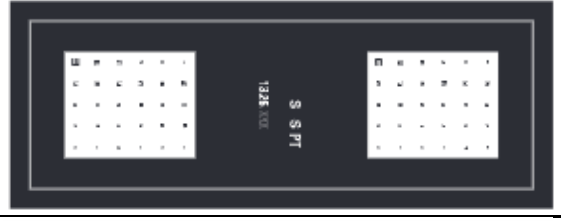
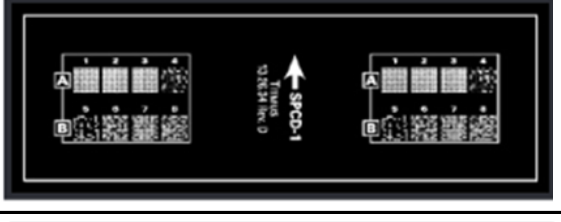

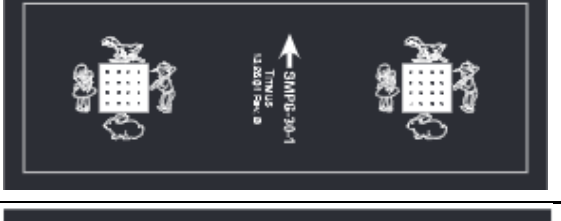

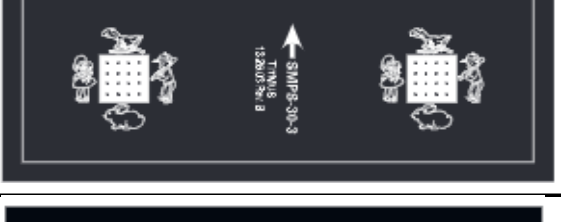


Titmus V3 Illinois Vision Screener

Slide Package (SH # 52784)

Slide	Slide Name	Testing Function	Testing Function Detail	Image
1	SPS-3 Muscle Imbalance	Muscle Balance	This is a Muscle Balance Test to test Vertical and Lateral Phoria in children. For children that fail this test, it helps identify if there is a possibility of amblyopia (lazy eye) or strabismus (squint eyes). Test is administered first at far and then at near.	
2	SAS-PT "E" Symbol Acuity	Acuity, Right & Left Eye	Tests visual acuity (sharpness of the eye). Right and Left Eye images are different. Administer test at FAR and using PLUS lenses. Refer to page 6 for Plus Lens Test.	
3	SPCD-1 Pediactric Color	Color Perception	Tests color perception. Test uses photographically reproduced TITMUS Pseudo-Isochromatic blocks.	
4	SA-BRL-4 Visual Acuity	Acuity, Both, Right & Left Eye	This test checks for visual acuity of Both Eyes, Right Eye and Left Eye using the same slide.	
5	SMPS-30-1 MI Preschool Acuity	Acuity, Right & Left Eye	Visual acuity testing using the Michigan Preschool Test. Right and Left Eye Images are different. Administer test at FAR. This is a part of a set of 4 slides.	
6	SMPS-30-2 MI Preschool Acuity	Acuity, Right & Left Eye	Visual acuity testing using the Michigan Preschool Test. Right and Left Eye Images are different. Administer test at FAR. This is a part of a set of 4 slides.	
7	SMPS-30-3 MI Preschool Acuity	Acuity, Right & Left Eye	Visual acuity testing using the Michigan Preschool Test. Right and Left Eye Images are different. Administer test at FAR. This is a part of a set of 4 slides.	
8	SMPS-30-4 MI Preschool Acuity	Acuity, Right & Left Eye	Visual acuity testing using the Michigan Preschool Test. Right and Left Eye Images are different. Administer test at FAR. This is a part of a set of 4 slides.	